

NEW

Micro Mixer E-022T

This well plate/tube mixer supports routine ELISA work and various microtubes with agitation speed of 300–2500 r/min and quiet operating sound.

•High-end model with temp. control --> P.022





Features

- Optimum for Stirring of 96-well plates and Microtubes
- The rack for Microtubes is available as an option


Applications

- Promotion of Eenzyme digestion of mouse tail, etc.
- Promotion of ELISA reaction, Biopanning, etc.
- Plasmid DNA extraction

Outline drawing and dimensions of Well and Deep Well

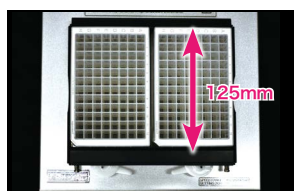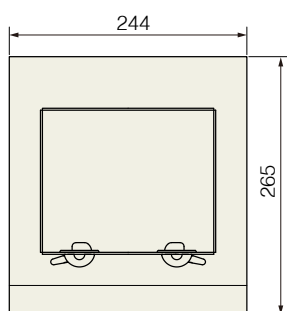

For the E-022T and the MBR-022R/K, if the long side of the Well is wider than 125 mm, use the fixing adapter that is supplied with the main unit.

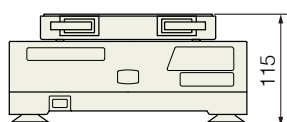

Microtube rack and Shaking rack for the E-022

Stacked frame exclusively for Next series-type. Countermeasures against earthquake ground motion have made it possible to install them in two stages.

It can be used for saving space in the laboratory, such as space on the equipment and effective use of dead space under the desk. Also, please use it in the temperature chamber. (Three stacks can be stacked for the E-022 as shown on the previous page.)

※For dimensions, refer to the next page.



SKE-2

SKE-3

SKT-2

(*) E-022T cannot be installed on the top row.

Model	E-022T
Shaking motion	Horizontal eccentric vibration (S-Max drive)
Vibration speed	300 to 2500 r/min
Shaking width	2mm
Speed display/Setting	Digital display/Input
Ambient temp. range	+4°C to +40°C
Applicable vessels	Well plate (2 pcs. Deep well can be installed as standard 2 pcs.) Microtube (requires P.069 optional racks)
Safety devices/ functions	Fuse, Overvoltage/overcurrent protection, Nonvolatile memory error, Speed limit error, Display/measuring range error, Motor overload
Dimensions (W × D × H)	244×265×115 Hmm
Weight	Approx. 6 kg
Standard accessories	Fuse (with built-in Fuse Holder) × 1 pc, Plate adapter
Power supply	AC 100 to 240 V

Product Name/Model	Capacity of vessels/Remarks
SKE-2	The E-022 can be loaded up to 2 units. 361×330×205 Hmm, Weight: Approx. 3.0 kg
SKE-3	The E-022 can be loaded up to 3 units. 361×439×393 Hmm, Weight: Approx. 5.6 kg
SKT-2	The E-022T can be loaded up to 2 units. 307×406×235 Hmm, Weight: Approx. 約3.5 kg

Optional: Microtube rack



EM-0524



EM-1524



EM-1515

Product Name/Model	Capacity of vessels/Remarks
Microtube rack EM-0524	0.5 mL×24 Dimension:126×86×28.5 Hmm, Weight : 80 g
Microtube rack EM-1524 (*1)	1.5/2.0 mL×24 Dimension:126×86×40.5 Hmm, Weight : 105 g
Microtube rack EM-1515 (*1,2)	5.0 mL×12 Dimension:126×86×55.5 Hmm

(*1) Watson-type 5 mL tubing is suitable. Eppendorf type does not fit.

(*2) Use the shaking rate at 1500 r/min or less.